

Migration analysis of perfluorinated compounds (PFC) in real food samples

Ludwig Gruber¹, Martin Schlummer¹, Jan Ungewiss¹, Gerd Wolz¹, Angela Möller¹, Nicole Weise¹, Hermann Fromme²

¹Fraunhofer Institute for Process Engineering and Packaging (IVV), Giggenhauser Straße 35, 85354 Freising, Germany, email: gruber@ivv.fraunhofer.de, phone: ++49 8161 491 706

²Department of Environmental Health, Bavarian Health and Food Safety Authority, Veterinärstraße 2, D-85764 Oberschleißheim, Germany

Introduction

There are indications, that food may play a major role as source for perfluorinated compounds (PFC) for the non-occupationally exposed population. One possible source for PFC in food is the migration from food contact materials. Begley and coworkers [1] demonstrated that the traditional food simulants water and oil are not the both extremes for these tensidic compounds. Furthermore the presence of soy emulsifier (lecithin, 0.05%) can change oil into a potent solvent for migration of PFC. Therefore it is recommendable to perform these migration analyses not in food simulants but in real food samples. The aim of this study was to develop a robust analytical method for different kinds of food samples. The development of the method was based on a publication by Taniyasu et al.

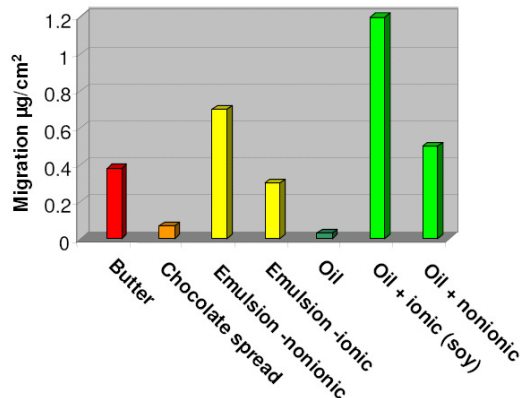


Figure 2: Example for PFOA migration from coated paper [3]

Method

Sample preparation:

The food sample has to be blended and homogenized. A 30 g aliquot of the homogenized sample is filled in a PP tube and fortified with internal standards. The internal standards ¹³C₄-PFOA and ¹³C₄-PFOS were purchased from Wellington (Ontario, Canada).

After addition of 2 ml of water and 10 minutes of ultrasonic extraction, the samples are centrifuged for 10 minutes at least with 4500 rpm and the supernatant has to be collected in a PP vessel. The sediment has to be twice re-extracted with methanol in an ultrasonic bath.

The combined extracts are diluted with the fivefold amount of water, filtrated through a syringe filter and subjected to solid phase extraction (SPE). Weak anion exchange SPE cartridges (Oasis® Wax, 150 mg, 6 ml, 30 µm) are preconditioned with 2 ml of methanol and water, respectively, and the extracts are passed through the preconditioned cartridges. The cartridges are then washed with methanol/water (1/1; v/v) and eluted with 1% NH₄OH in methanol. SPE eluates were evaporated under a gentle stream of nitrogen and diluted with water to a final volume of 1 ml.

HPLC-ESI-MS/MS:

Identification and quantification of perfluorinated substances was performed on a Surveyor Plus HPLC connected to a Quantum Ultra AM mass spectrometer (both Thermo Fisher). Chromatographic separation was achieved by a Fusion RP phase (20 x 2 mm, 2 µm, Phenomenex). Gradient HPLC was performed with methanol and 5mM ammonia acetate in water (pH 3.5), increasing methanol from 20 to 100% within 10 minutes. Mass spectrometry was performed by electron spray ionization (ESI) in the negative ion mode and subsequent single reaction monitoring.

Even with HPLC-ESI-MS/MS interferences can disturb the quantification of PFOS significantly [3]. Taurodeoxycholate, which coelutes with PFOS isomers, produces significant interference on the most sensitive MS/MS transition of PFOS, 499 to 80. The interference can be avoided by monitoring the more specific transition 499 to 99 [4].

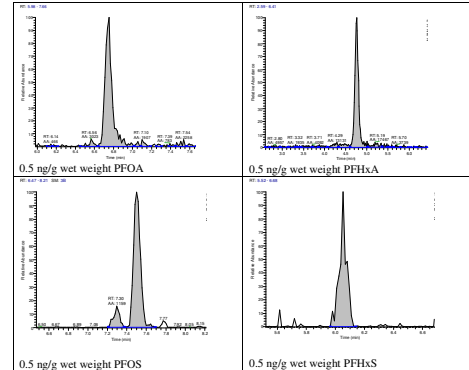


Figure 2: Sensitivity of the method shown by MS/MS fragmentograms of PFOA, PFOS, PFHxA, and PFHxS reflecting 0.5ng/g wet weight in diet samples

Results and Conclusions

The described method proved to be a robust and sensitive means to analyze PFC in food. In fish and mixed food samples, limits of detection were between 0.05 ppb for PFOS and PFOA and 0.1 ppb for the other target compounds. This analytical method enables laboratories to perform migration experiments for perfluorinated compounds (PFC) with real food samples in order to avoid problems for exposure estimates.

References

- [1] Begley TH, White K, Honigfort P, Twaroski ML, Neches R, Walker RA (2005). Perfluorochemicals: potential sources of and migration from food packaging. *Food Additives and Contaminants* 2005; 22:1023–1031.
- [2] Taniyasu S, Kannan K, So MK, Gulkowska A, Sinclair E, Okazawa T, Yamashita N (2005). Analysis of fluorotelomer alcohols, fluorotelomer acids, and short- and long-chain perfluorinated acids in water and biota. *Journal of Chromatography A* 2005; 1093:89
- [3] Begley TH, Hsu W, Noonan G, Diachenko G (2008). Migration of fluorochemical paper additives from food-contact paper into foods and food simulants. *Food Additives and Contaminants* 2008; 25: 384–390.
- [4] Gruber L, Schlummer M, Ungewiss J, Wolz G, Möller A, Weise N, Fromme H (2007). Analysis of sub-ppb levels of perfluorooctanoic acid (PFOA) and perfluorooctanesulfonate (PFOS) in food and fish. *Dioxin 2007, 27th International Symposium on Halogenated Persistent Organic Pollutants : 2.-7. September 2007, Tokyo, Japan, (Organohalogen compounds 69)*
- [5] Fromme H, Schlummer M, Möller A, Gruber L, Wolz G, Ungewiss J, Böhmer S, Dekant W, Mayer R, Liebl B, Twardella D (2007). Exposure of an adult population to perfluorinated substances using duplicate diet portions and biomonitoring data. *Environmental Science and Technology* 2007; 41:7928-7933

Acknowledgements

This study was partly financed by the Bavarian State Ministry of the Environment, Public Health and Consumer Protection for financial support. The authors thank all members of the INES study group and especially the colleagues of the Department of Environmental Health of the Bavarian Health and Food Safety Authority and the Fraunhofer IVV for excellent cooperation.



Fraunhofer Institut
Verfahrenstechnik
und Verpackung